



WHITE · PEARL

BY OSSEUS

PRECISE PLACEMENT

Patented cannulated design allows for quick and accurate preliminary fixation at midline of plate.

EXTREME ANGULATION


System framework offering $\pm 20^\circ$ cephalad / caudal angulation allows for maximum flexibility on challenging alignments.

ROBUST LOCKING

Manual, one-step locking rivets provide a visual, tactile, and audible confirmation that screws have been locked.

LOW PROFILE


Interchangeable Ultra-thin plate designed for decreased anatomic distribution and increased disk space visibility.

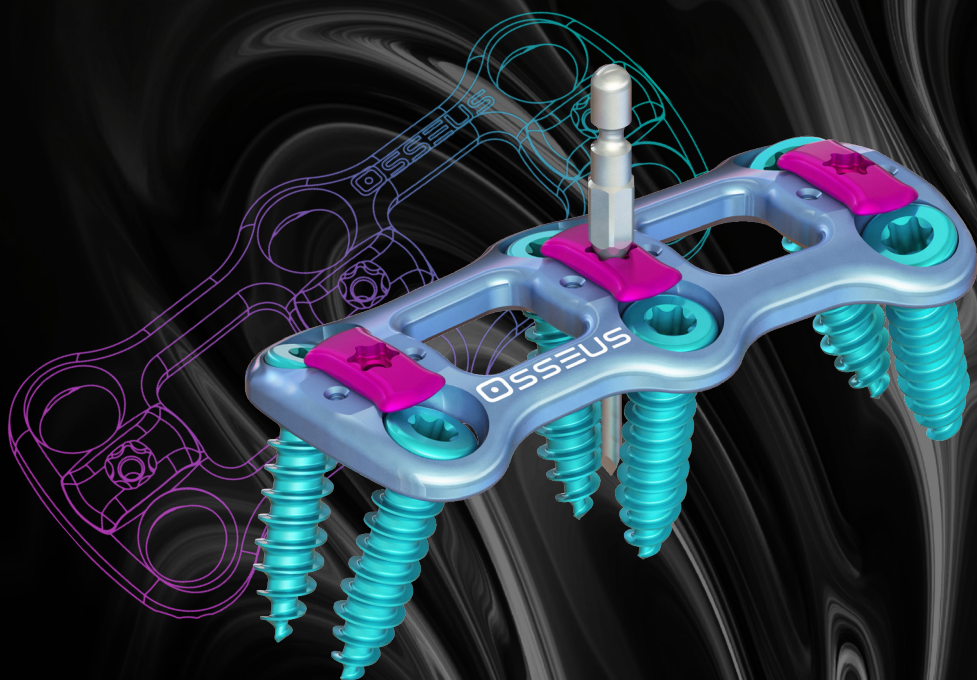


The White Pearl ACP is a pioneering cervical spine fixation system that is ideal for the treatment of degenerative or traumatic cervical disorders. The plate's design allows for quick and accurate preliminary fixation as well as maximum screw angulation to accommodate challenging anatomy. Numerous screw options accommodate all preferences and are backed up by robust manual locking rivets to prevent screw expulsion.

TECHNICAL SPECIFICATIONS

WHITE PEARL CERVICAL PLATES

	1 LEVEL	2 LEVEL	3 LEVEL	4 LEVEL
	8.5mm	24mm	39mm	60mm
	10mm	26mm	42mm	64mm
	12mm	28mm	45mm	68mm
	14mm	30mm	48mm	72mm
	16mm	32mm	51mm	76mm
	18mm	34mm	54mm	80mm
	20mm	37mm	57mm	84mm
	22mm	40mm	60mm	
	24mm	43mm	63mm	
	26mm	46mm	66mm	
			69mm	

WHITE PEARL BONE SCREWS

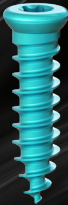	STYLES:	DIAMETERS:	LENGTHS:
	Fixed	Ø3.75mm	10mm
	Variable	Ø4.25mm	12mm
	Self - Drilling		14mm
	Self - Tapping		16mm
			18mm
			20mm

WHITE PEARL PRELIMINARY FIXATION

PIN	SCREW	TACK
		